




be free Day Pro - Product Specifications



Be Free Day Pro is an individually crafted silicone hydrogel contact lens specifically designed for Myopia Management. It is powered by the Brien Holden Vision Institute's patented Extended Depth of Focus (EDOF) technology, which slows myopia progression and supports a comfortable adaptation to the lens, enhancing the overall wearing experience. A monthly disposable contact lens, BeFreeDay features high water content, low coefficient of friction and low elastic modulus, which combine to improve comfort throughout the day. Its wide range of parameters ensure an excellent fit, especially for the youngest contact lens wearers.

PARAMETERS	
Base curves (mm)	7.10 to 9.80 (0.30 steps)
Diameters (mm)	13.50 to 15.50 (0.50 steps)
Spheres (D)	-0.25 to -15.00 (0.25 steps)
Cylinders (D)	-0.75 to -8.00 (0.25 steps)
Axes (°)	All (1° steps)

MATERIAL	
Type	Filcon 5b (60) [75%]
DK (ISO 9913-1-1998)	60
DK/t (-3.00 D)	50
Water Content	75%
Central Thickness (-3.00 D)	0.12
Cof	0.02
Modulus	0.33
UV filter	Class 1
Handling tint	Blue
Pack size	3 Lenses
Manufacturing Process	Lathed

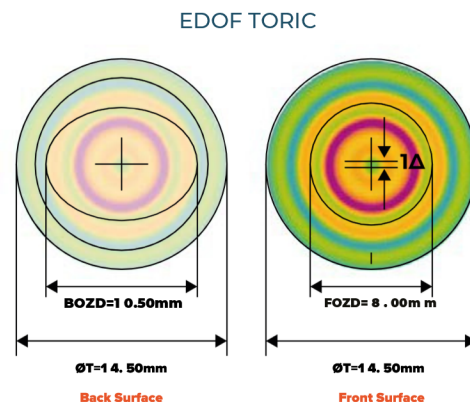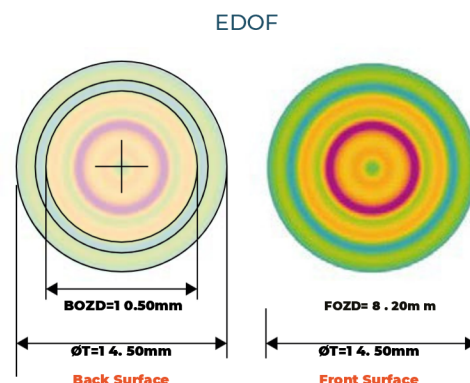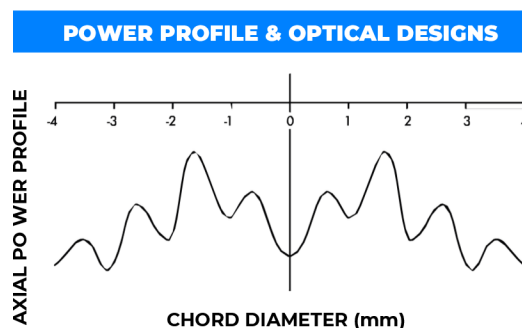

Find out more about
Be Free Day Pro
by visiting our product page.

